

Physics 1250

Laboratory Activities & Worksheets

Sixth Edition

Department of Physics
The Ohio State University

Copyright 2021 by the Department of Physics, The Ohio State University

All rights reserved.

Permission in writing must be obtained from the publisher before any part of this work may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying and recording, or by any information storage or retrieval system.

Printed in the United States of America

10 9 8 7 6 5 4 3 2 1

ISBN 978-0-7380-7642-3

Hayden-McNeil Publishing
14903 Pilot Drive
Plymouth, MI 48170
www.hmpublishing.com


TABLE OF CONTENTS

Lab Instructions

Introduction to *LoggerPro*

Experiment I -	1-D Kinematics
Experiment II -	Vectors
Experiment III -	2-D Kinematics
Experiment IV -	Dynamic Forces
Experiment V -	Static Friction
Experiment VI -	Conservation of Energy
Experiment VII -	Conservation of Momentum
Experiment VIII -	Energy and Momentum
Experiment IX -	Rotational Dynamics
Experiment X -	Vibrations
Experiment XI -	Fluids
Experiment XII:A -	Specific Heat
Experiment XII:B -	Ideal Gas
Experiment XII:C -	Heat Engine
Experiment XIII -	Special Relativity

Worksheets

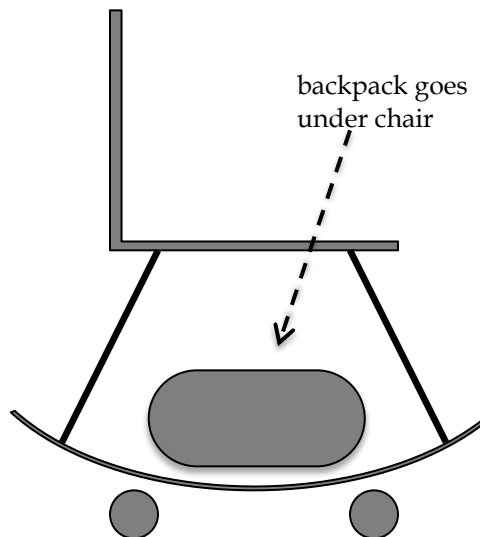

Lab Instructions

The two-hour lab session lies at the heart of this course. You will work a combination of problems and laboratory activities in small groups with the assistance of the instructor. These sessions are where you can best learn the material, in order to be able to complete the homework assignments and be prepared for the quizzes and exams that will determine your grade.

Typically, you and your group will be given a combination of laboratory activities and problems to do in a set order. Most activities and problems are not graded for accuracy, though some of the laboratory activities may be graded (and are so marked). Your “lab score” is based primarily on participation with your group and learning. You are done when either your lab instructor determines that you have finished or the lab time has elapsed.

It is better that you work a few problems well with complete understanding than rush through a large number without understanding your work. You may not always finish all of the problems.


Your lab room may be equipped with whiteboards or blackboards. You are encouraged, and may be required, to use them for the group problems and/or for lab activities.

You should record all measurements and calculations in your individual lab manual (all members of your group are expected to do so). Your instructor may ask to examine it at any time. You are encouraged to ask for help from your instructor at any time. You are welcome to use your textbook and class notes. Also, do not forget that you are part of a group and are expected to work and learn together. Bring this manual and a scientific calculator to class.

Be on time. If you come late, you may not be allowed to perform the lab.

Academic Misconduct: (These considerations are in addition to the general University rules.) All laboratory work will be done during the laboratory period. Do not bring completed lab work to the laboratory. The presence of completed material in the lab book can be considered evidence of intent to commit academic misconduct. Do “graded activities” only when your instructor is watching. Do not mistreat the laboratory equipment -- in particular, do not write on the apparatus. If you are uncertain what is permitted while doing a given graded lab activity, ask your instructor before proceeding.

Introduction to LoggerPro

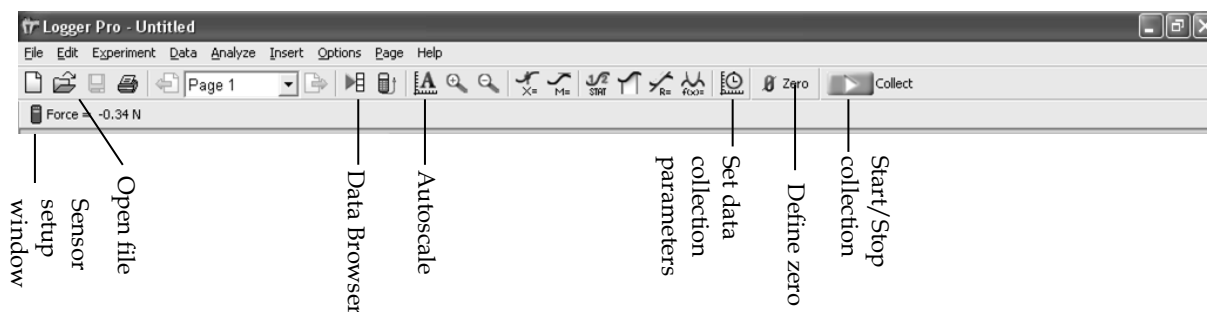

To start data collection, click the green “Collect” button on the tool bar. There is a delay of a second or two before data collection starts, so don’t try to time clicking the button with your actions.

Many parameters in the software are adjustable, like the sample collection time. Adjusting them is done through the menus, but some common commands and parameters are available through the toolbar. For example, click the clock icon to change the sample rate or time.

Often you will be told to load/open a file. Such files will set all the relevant parameters for the experiment.

Graph parameters can be changed by clicking on the graph to select it and using the menus, or by right-clicking on the graph. Common graph commands include *zoom* and *autoscale*.

If you can't see your data (the graphs appear blank) or the data is a flat line at zero you probably need to adjust the graph scaling. Start with “Autoscale from 0”, then zoom-in or adjust the axes.

If the “collect” button is gray and can't be clicked then there may be a problem with the sensor hardware. Make sure the cable to the sensor is snapped in place and likewise with the cables going to the interface box. (The cable from the sensor goes to the interface box; just follow it.) If the button is still gray, the problem might be the software. Quit out of LoggerPro and restart it.

If all else fails, ask your lab instructor for help.